**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS – DO NOT DESTROY	3	7. Washing instructions.....	8
2. Introduction	6	8. Storage and Disposal.....	8
3. Package Contents	6	9. Troubleshooting Guide.....	9
4. Parts and Controls	6	10. FCC Compliance Information	9
5. Initial Use	7	11. Product Specifications.....	10
6. To Use	7	12. Warranty	10

IMPORTANT INSTRUCTIONS - DO NOT DESTROY SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service at 1-800-536-0366 or at info@beurer.com.

Dear Valued Customer

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer Team

1. IMPORTANT SAFETY INSTRUCTIONS – DO NOT DESTROY

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

WARNING

Failure to follow instructions may cause overheating, fire or personal injury.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 WARNING	WARNING: Indicates a hazardous situation which, if not avoided, could result in death or serious injury.
 CAUTION	CAUTION: Indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
 NOTICE	NOTICE: addresses practices not related to personal injury, such as product and/or property damage.

WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- Be sure to read carefully the accompanying instruction booklet before using this blanket.
- This blanket is only designed to warm up the human body. Do not use this product for any other purpose.
- Keep children away from packaging materials (risk of suffocation).

- Do not use this blanket with an infant, an immobile person, or anyone insensitive to heat, such as a person with poor blood circulation.
- Never leave the blanket unattended
- Never pull this blanket by the power cord. Do not use the cord as a handle. Unplug when not in use.
- Do not use the blanket when symptoms of appendicitis are present. Consult your physician.
- Do not use an electrical outlet that has become loose or does not engage the power plug tightly.
- Do not remove product labels that contain warnings or safety instructions.
- Consult with a physician for use with controlled substances.
- Do not use on areas of insensitive skin. This may lead to skin burns.
- Burns can occur regardless of the control setting. Check skin under the blanket frequently.
- Do not use the blanket with liniment, salve or ointment preparations that contain heat-producing ingredients. Skin burns can result.
- Do not use with alcohol.
- Avoid sharp folds.
- Do not tuck in the wired area of this blanket.
- Do not route electrical cords between the mattress and box spring. Do not allow the cords to be pinched.
- Keep the controller away from drafty areas that may become damp or wet such as an open window.
- Do not use this blanket if you or the product is wet.
- Do not use pins – They may damage the electric wiring.
- Be sure that this blanket is used on an alternating-current supply circuit of the proper voltage.
- Turn off or unplug when not in use.
- Excessive heating may result from folding blanket, other than a mattress, when in actual use.

- IF IMPROPER OPERATION OF THIS BLANKET IS OBSERVED, DISCONTINUE ITS USE IMMEDIATELY AND CONSULT THE MANUFACTURER OR THE DEALER REGARDING ITS REPAIR.
- Do not use with rocking chairs, recliners, chairs with casters or wheels, or other furniture subject to movement.
- Do not tuck the blanket or supply cord in a sofa or other furniture. All portions of the blanket or controller should be visible at all times.
- Wrapping of cord around the controller may damage the cord. Loop loosely when storing.
- Do not trap or pinch the cord.
- NOT FOR HOTEL, MOTEL OR OTHER INSTITUTIONAL USE.
- This blanket has a polarized plug (one blade is wider than the other) as a safety feature. This plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.
- Do not dry-clean or bleach this blanket. Cleaning or bleaching solvents may have a deteriorating effect on the insulation of the heating element.

 **CAUTION – To reduce the risk of personal injury or product/property damage:**

- This product is intended solely for household use and not for medical or commercial purposes.
- Check this blanket and cord frequently for signs of wear and tear or damage. If detected, discontinue use and contact customer service.

 **NOTICE**

- The device should be cleaned from time to time. Do not use any abrasive cleaning products and never submerge the blanket in water.
- Do not place any objects on the blanket when not in use.

- Protect the blanket from knocks, damp, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).

2. Introduction

We designed this blanket to keep you warm while lounging in the comfort of your home. It's also great to take with you when traveling.

3. Package Contents

- 1 x HD 71 Heated Blanket
- 1 x Controller
- 1 x Instruction manual

4. Parts and Controls


Safety System

NOTICE

The blanket is equipped with a safety system. The “high” indicator light will blink red when an error occurs. If this happens, unplug the blanket from the electrical outlet. Check the physical condition of the cord, the wall outlet, and the connection between the Control Connector and the Blanket Connector. If blinking persists after you confirm a proper connection, email info@beurer.com.

5. Initial Use

Attach the Blanket Connector to the Control Connector. Push both parts together securely to ensure a tight fit.

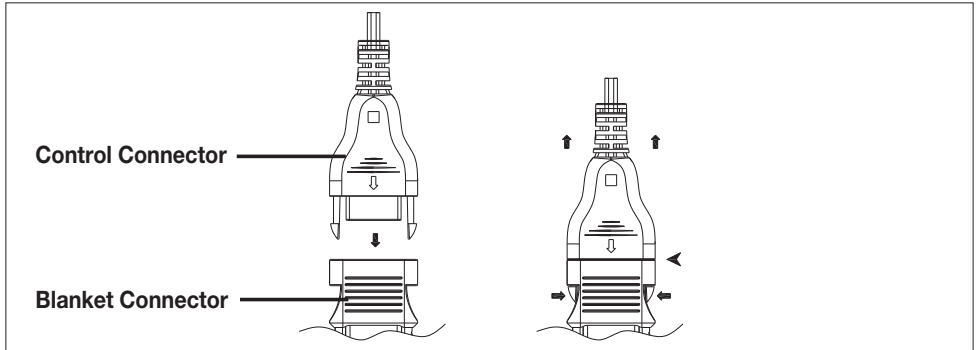i NOTICE

The Control Connector must be connected to the Blanket Connector before plugging the unit into an electrical outlet. After connecting the Blanket Connector to the Control Connector, plug the unit into a 120V AC electrical outlet.

6. To Use

To turn the blanket on, press the power button. Select the desired heat setting by continuing to press the power button until the heat indicator light is the color of your desired heat setting:

Green=Warm Orange=Medium
Yellow=Low Red=High

To turn the blanket off, press the power button one more time.

i NOTICE

If the unit does not turn on when you press the power button, make sure the Control Connector is securely plugged into the Blanket Connector.

To turn the blanket off, scroll through all the heat settings until the heat indicator light turns off.

! WARNING

If the blanket is being used over several hours, we recommend that you set the lowest temperature setting on the controller to avoid overheating the heated body, which may lead to burns to the skin.

AUTO-OFF FEATURE

To conserve energy and prevent skin burns caused by excessive use, the blanket includes an auto-off feature. The unit will turn off automatically after two hours of continuous use. The red auto-off light will blink when this feature is activated.

i NOTICE

After a power outage, the blanket automatically resets the timer on the auto-off feature.

7. Washing instructions

WARNING

Always remove the plug from the socket and Control Connector from Blanket Connector before cleaning the blanket. Otherwise there is a risk of electric shock.

Use a dry, lint-free cloth for controller cleaning and care.

CAUTION

Never immerse the controller in water or other liquids. It may otherwise be damaged.

CAUTION

Do not use chemical cleaning or scouring agents to clean the blanket and controller. The blanket and controller may otherwise be damaged.

CAUTION

Risk of damage to the product or property! Do not dry clean! Do not use bleach! Do not use Wringer!
Do not iron your blanket!

To prolong the life of your heated blanket, proper care and cleaning is essential. Follow these washing instructions to preserve the softness of the fabric and protect the electrical components:

Wash by hand or as follows in automatic washer using slow agitation and normal spin speed: Fill washer with lukewarm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure blanket before washing. Add blanket and soak five minutes. Agitate two minutes. Spin. Fill for cool rinse, agitate one minute then spin. Dry in automatic dryer on low heat, two minutes only, or line dry. Remove damp blanket from dryer, block to original size, and drape to finish drying. Do not iron your blanket.

CAUTION

Make sure the Blanket Connector is completely dry before plugging the Blanket into an electrical outlet.

8. Storage and Disposal

- Allow the heating seat to cool down before storing.
- To prevent product damage do not place any objects on it.
- If you do not plan to use the heating seat for a long period of time, store it in the original packaging in a cool, dry place.
- Disconnect the Control Connector and the controller for storage.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

9. Troubleshooting Guide

Problem	Cause	Solution
<p>The indicator lamps are not illuminated while</p> <ul style="list-style-type: none"> • the controller is properly connected to the blanket • the mains plug is connected to a working socket • the power button for ON/OFF and temperature settings is set to “Warm”, “Low”, “Medium” or “High” 	<p>The safety system has switched off the blanket permanently.</p>	<p>Contact Customer Service for assistance.</p>
<p>The “High” indicator light blink red.</p>	<p>An error occurs.</p>	<p>Unplug the blanket from the electrical outlet. Check the physical condition of the cord, the wall outlet, and the connection between the Control Connector and the Blanket Connector. If blinking persists after you confirm a proper connection, contact Customer Service for assistance.</p>

10. FCC Compliance Information

HD 71 Heated Blanket

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: 1. This device may not cause harmful interference, and 2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user’s authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. Product Specifications

Blanket Size:	50 x 60 in (130 x 150 cm)
Power Supply:	AC: 120V, 60Hz, 100W
Blanket Material:	Microplush
Auto Shut-Off Timer:	120 minutes
4 Heat Settings:	Warm-temperature range Low-temperature range Medium-temperature range High-temperature range

Subject to technical modifications.

12. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Heated Blanket, model HD 71, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners. We will, at our option, repair or replace the Beurer Heated Blanket, model HD 71, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement. For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty. **THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you. **IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you. Beurer does not authorize anyone, including, but not limited to, retailers, the

subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein. This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers.

Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer. This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

Made in China

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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