|                 |                                     |    |
|-----------------|-------------------------------------|----|
| <b>ENGLISH</b>  | <b>Handheld Massager</b>            |    |
|                 | Instructions for use .....          | 3  |
| <b>ESPAÑOL</b>  | <b>Masajeador Portátil</b>          |    |
|                 | Manual de instrucciones .....       | 8  |
| <b>FRANÇAIS</b> | <b>Appareil de Massage Portatif</b> |    |
|                 | Mode d'emploi.....                  | 14 |

Distributed by / Distribuido por /  
Distribué par:  
Beurer North America LP  
1 Oakwood Blvd, Suite 255  
Hollywood, FL 33020, USA  
[www.beurer.com](http://www.beurer.com)

-  **Questions or comments?**  
Call our US-based customer service toll free at 1-800-536-0366.
-  **¿Preguntas o comentarios? Llame a nuestro servicio al cliente en Estados Unidos gratis al teléfono 1-800-536-0366.**
-  **Questions ou commentaires? Appelez gratuitement notre service clientèle situé aux États-Unis au 1-800-536-0366.**

**READ THIS MANUAL COMPLETELY AND CAREFULLY  
BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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## CONTENT

|                                       |   |                                         |   |
|---------------------------------------|---|-----------------------------------------|---|
| 1. IMPORTANT SAFETY INSTRUCTIONS..... | 3 | 5. To Use .....                         | 6 |
| 2. Introduction .....                 | 5 | 6. Care, Maintenance and Disposal ..... | 6 |
| 3. Package Contents .....             | 5 | 7. Technical Specification .....        | 6 |
| 4. Parts and Controls.....            | 6 | 8. Warranty .....                       | 7 |

## SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

### Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support

Visit [www.registerbeurer.com](http://www.registerbeurer.com) or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With best regards,

Your Beurer Team

## 1. IMPORTANT SAFETY INSTRUCTIONS

### Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

**READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.**

|                                                                                                   |                                                                                                                                                                          |
|---------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|                 | This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death. |
|  <b>DANGER</b>  | <b>DANGER</b> indicates a hazardous situation which, if not avoided, will result in death or serious injury.                                                             |
|  <b>WARNING</b> | <b>WARNING</b> indicates a hazardous situation which, if not avoided, could result in death or serious injury.                                                           |

|                                                                                                |                                                                                                                 |
|------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------|
|  <b>CAUTION</b> | <b>CAUTION</b> indicates a hazardous situation which, if not avoided, could result in minor or moderate injury. |
| <b>NOTICE</b>                                                                                  | <b>NOTICE</b> addresses practices not related to personal injury, such as product and/or property damage.       |
|  <b>NOTE</b>   | <b>NOTE</b> provides important information.                                                                     |

**DANGER** – To reduce the risk of electric shock:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Do not reach for an appliance that had fallen into water. Unplug immediately.
- Do not use while bathing or in a shower.
- Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

**WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons:

- For household use only.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it had a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return appliance to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the off position, then remove plug from outlet.
- The device has a heated surface. Persons sensitive to heat must be careful when using the device.
- Do not use massager in close proximity to loose clothing or jewelry.
- Keep long hair away from massager while in use.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Do not use this massager:
  - On children.
  - On animals.
  - On persons suffering from a herniated disk or from abnormal changes to the vertebrae.
  - To massage the heart region of persons fitted with a pacemaker. Always consult a doctor before using this massager on a person fitted with a pacemaker.
  - On swollen, burned, inflamed or damaged skin or body parts.
  - During pregnancy.
  - On the face/eyes, larynx, or other particularly sensitive parts of the body.
  - While asleep.
  - After taking medication or consuming alcohol.

- This device is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking an understanding of the equipment and how it is used.
- This device is not a toy; keep away from children.
- Consult your doctor before using the massage device if you suffer from:
  - A serious illness or are recovering from an operation.
  - Common cold accompanied by a fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica) or acute inflammations.
- Do not operate the device:
  - With a voltage different than that specified on the device.
  - If it or any of its accessories is damaged.
  - During a thunderstorm.
- In the event of damage or a malfunction, switch the device off immediately and disconnect it from the power supply.
- Never pull on the power cable or on the appliance in order to disconnect the plug.
- Never hold or carry the device by the power cable. Keep the cables away from hot surfaces.
- Never open the device for any reason.
- If the power cord of the device is damaged, have it replaced right away.
- Make sure the massage device, plug and the power cord do not come into contact with water or other fluids.
- Operate the device:
  - Only in dry inside rooms
  - Only with dry hands
  - Never in the bath, in the shower or near a filled basin
  - Never in swimming pools, whirlpools or in the sauna
  - Never out in the open
- Never try to retrieve the device if it falls into water. Disconnect the plug immediately from the electrical outlet.
- Repairs to electrical appliances must only be carried out by qualified persons.
- Never leave the device unsupervised, particularly if children are around.
- Never use the device under covers, such as blankets or pillows.

#### **CAUTION:**

- Switch the device off and disconnect the power supply after every use and before every cleaning.
- Make sure that all moving parts can move freely at all times.
- Make sure hair cannot get caught in the device.
- Avoid contact with sharp or pointed objects.
- Do not kink the power cord.
- Continuous operation for no more than 15 minutes.
- Make sure the device is not too close to heat sources such as ovens, radiators or direct sunlight.
- To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not change the plug in any way.


## **2. Introduction**

This electric handheld massager allows you to give yourself or other people soothing massages. Massages can be relaxing or stimulating. This massager delivers powerful and intensive back massages and is intended for private use only and not for medical or commercial purposes.

## **3. Package Contents**

1x MG 21 Handheld Massager  
 3x Interchangeable massage attachments  
 1x Instruction Manual

## 4. Parts and Controls
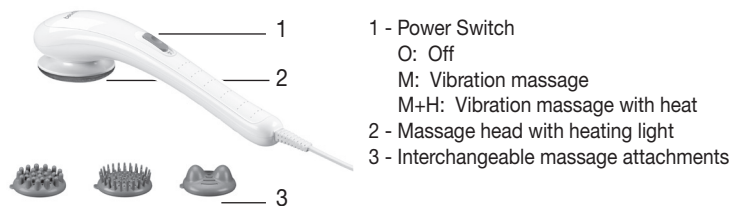
## 5. To Use

Attach one of the supplied Massage Attachments over the massage head. Make sure the device is switched off and then plug the power cord into a standard electrical outlet. Switch the device on by sliding the power switch to the desired massage setting “M” or “M+H.”

Guide the device with both hands to the area of the body to be massaged. You can stand, sit, or lie down. You can alter the massage intensity by pressing down firmly or gently. Do not massage one part of the body longer than three minutes before moving to other areas.

### **WARNING:**

- Stop the massage or change the position of the massager if it is painful or unpleasant.
- Do not jam or force any part of your body, particularly your fingers, into the massage head.

### **NOTE**

Do not use the massage device right before going to bed. The massage can have a stimulating effect and may make it difficult to fall asleep.

## 6. Care, Maintenance and Disposal

### **WARNING:**

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Ensure no fluids enter the device or accessories.
- Do not use the device until it is completely dry.

### **Cleaning**

Unplug the device and allow it to cool before cleaning it with a damp cloth and, if necessary, a mild detergent.

### **Storage**

Keep the massage device in a cool and dry place where it is out of the reach of children, preferably in its original packaging.

### **Disposal**

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

## 7. Technical Specification

|               |                                                |
|---------------|------------------------------------------------|
| Type          | MG 21                                          |
| Dimensions    | 2.5 in W x 3.1 in H x 13 in (65 x 80 x 330 mm) |
| Weight        | 0.9 lbs (0.41 kg)                              |
| Rated output  | 7 W                                            |
| Rated voltage | 120 V~ 60 Hz                                   |

## 8. Warranty

### Limited Lifetime Warranty For Original Purchaser

Your Beurer handheld massage, model MG 21, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer handheld massage, model MG 21, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

**THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: [www.beurer.com](http://www.beurer.com)

Distributed by:  
Beurer North America LP  
1 Oakwood Blvd, Suite 255  
Hollywood, FL 33020, USA  
[www.beurer.com](http://www.beurer.com)

**Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.**

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